

Human Properdin

HK334

Edition 10-13

**ELISA KIT
PRODUCT INFORMATION & MANUAL**

Read carefully prior to starting procedures!
For use in laboratory research only
Not for clinical or diagnostic use



Note that this user protocol is not lot-specific and is representative for the current specifications of this product. Please consult the vial label and the Certificate of Analysis for information on specific lots. Also note that shipping conditions may differ from storage conditions.

For research use only. Not for use in or on humans or animals or for diagnostics. It is the responsibility of the user to comply with all local/state and federal rules in the use of this product. Hycult Biotech is not responsible for any patent infringements that might result from the use or derivation of this product.

TABLE OF CONTENTS

	Page
1. Intended use	2
2. Introduction	2
3. Kit features	2
4. Protocol overview	3
5. Kit components and storage instructions	5
Materials requires but not provided.....	5
6. Warnings and precautions	6
7. Sample preparation	7
Collection and handling	7
Dilution procedures.....	7
8. Reagent preparation.....	8
9. ELISA protocol	10
10. Interpretation of results.....	11
11. Technical hints	11
12. Quality control	11
13. Performance characteristics	12
14. Troubleshooting	13
15. References.....	13

1. INTENDED USE

The human Properdin ELISA kit is to be used for the *in vitro* quantitative determination of human Properdin in serum, plasma, bronchoalveolar lavage fluid, urine and cerebrospinal fluid samples. This kit is intended for laboratory research use only and is not for use in diagnostic or therapeutic procedures.

The analysis should be performed by trained laboratory professionals.

2. INTRODUCTION

Properdin, or factor P, is a single-chain plasma glycoprotein (apparent molecular mass 52-55 kDa), consisting of six thrombospondin repeat sequences between short N- and C-terminal domains. In the blood it exists as a mixture of head-to-tail dimers, trimers and tetramers, with a preponderance of trimers.

The protein is produced by a variety of leukocytes, including monocytes, T lymphocytes and neutrophils, but also by endothelial cells in which properdin synthesis and release are induced by shear stress.

Properdin participates in the alternative pathway of complement activation together with C3 and factors B, D, I and H. It stabilizes the labile C3bBb which is deposited on immune complexes or foreign surfaces. This permits amplification of C3bBb formation in competition with catabolism of C3b by factor I, which uses factor H as a cofactor. The local amplification process leads to the creation of the alternative pathway C5 convertase, C3bBb3b, and initiates the terminal pathway of complement activation.

Properdin prolongs the half-life of surface-bound C3bBb from 1½ minute to about 18 minutes. This is due to several effects: inhibition of C3b cleavage by factor I, increased affinity for factor B and inhibition of the dissociation of C3bBb into C3b and Bb. Properdin is consumed by binding to C3bBb, this binding showing an order of preference of tetramers over trimers over dimers, which is also the order of functional activity of the oligomeric forms.

Deficiency or malfunction of the molecule may lead to severe impairment of alternative pathway activation, depending on the precise nature of the defect. One parameter of functional defect in the presence of measurable levels of protein is an impaired generation of the more active tetramer and trimer forms.

The properdin gene is located on the short arm of the X chromosome, and congenital properdin deficiencies are therefore inherited as typical X-linked recessive traits.

Three types of familial deficiency have been described: type 1 (or I) is characterized by serum with very low or absent properdin activity in hemolytic assays and <0.1 µg/ml immunoreactive protein; type 2 (or II) is characterized by low but detectable levels of immunoreactive protein (>2 µg/ml) and impairment of some, but not all functional test, and type 3 (or III) has normal levels of immunoreactive but dysfunctional protein. 10 Normal serum or plasma levels of properdin are variously quoted as about 5 or 25 µg/mL. Differences in properdin levels between healthy and diabetic persons have been determined (Chaudhry_B_2008).

3. KIT FEATURES

- Working time of 3½ hours.
- Minimum concentration which can be measured is 0.3 ng/ml.
- Measurable concentration range of 0.3 to 20 ng/ml.
- Working volume of 100 µl/well.

Cross reactivity


Potential cross-reacting detected in the human Properdin ELISA:

Cross reactant	Reactivity
Mouse	Average
Pig	Strong
Rabbit	Weak
Sheep	Negative
Horse	Negative
Goat	Negative
Rat	Negative

Table 1

Cross-reactivity for other species or proteins/peptides has not been tested.

4. PROTOCOL OVERVIEW

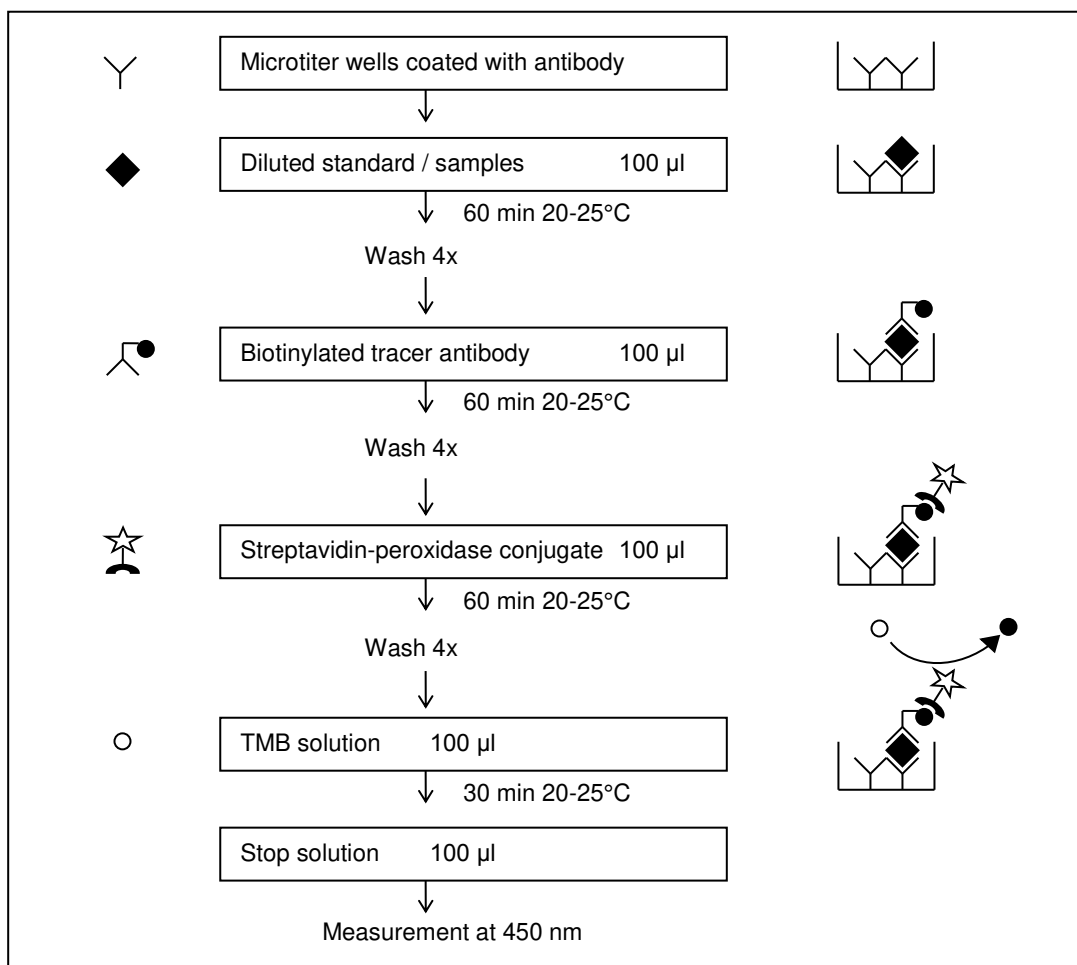

- The human Properdin ELISA is a ready-to-use solid-phase enzyme-linked immunosorbent assay based on the sandwich principle with a working time of 3½ hours.
- The efficient format of a plate with twelve disposable 8-well strips allows free choice of batch size for the assay.
- Samples and standards are incubated in microtiter wells coated with antibodies recognizing human Properdin.
- Biotinylated tracer antibody will bind to the captured human Properdin.
- Streptavidin-peroxidase conjugate will bind to the biotinylated tracer antibody.

- Streptavidin-peroxidase conjugate will react with the substrate, tetramethylbenzidine (TMB).
- The enzyme reaction is stopped by the addition of oxalic acid.
- The absorbance at 450 nm is measured with a spectrophotometer. A standard curve is obtained by plotting the absorbance (linear) versus the corresponding concentrations of the human Properdin standards (log).
- The human Properdin concentration of samples, which are run concurrently with the standards, can be determined from the standard curve.

5. KIT COMPONENTS AND STORAGE INSTRUCTIONS

Kit component	Cat. #	Quantity HK334-01	Quantity HK334-02	Color code
Wash buffer 20x	WB21	1 vial (60 ml)	1 vial (60 ml)	Colorless
Dilution buffer 10x	DB89	1 vial (20 ml)	1 vial (20 ml)	Green
Standard		2 vials, lyophilized	4 vials, lyophilized	White
Tracer, biotinylated		1 vial, 1 ml lyophilized	2 vials, 1 ml lyophilized	White
Streptavidin-peroxidase 100x	CON03	1 tube, 0.25 ml in solution	1 tube, 0.25 ml in solution	Brown
TMB substrate	TMB050/TMB100	1 vial, (11 ml)	1 vial, (22 ml)	Brown
Stop solution	STOP110	1 vial, (22 ml)	1 vial, (22 ml)	Red
12 Microtiter strips, pre-coated		1 plate	2 plates	
Certificate of Analysis		1	1	
Manual		1	1	
Data collection sheet		2	2	

Table 2

- Upon receipt, store individual components at 2 - 8°C. Do not freeze.
- Do not use components beyond the expiration date printed on the kit label.
- The standard and tracer in lyophilized form and the streptavidin-peroxidase in concentrated solution are stable until the expiration date indicated on the kit label, if stored at 2 - 8°C.
- The exact amount of the standard is indicated on the label of the vial and the Certificate of Analysis.
- The standard is single use. After reconstitution the standard cannot be stored.
- Once reconstituted the tracer is stable for 1 month if stored at 2 - 8°C.
- The streptavidin-peroxidase can only be stored in concentrated solution and cannot be stored once diluted.
- Upon receipt, foil pouch around the plate should be vacuum-sealed and unpunctured. Any irregularities to aforementioned conditions may influence plate performance in the assay.
- Return unused strips immediately to the foil pouch containing the desiccant pack and reseal along the entire edge of the zip-seal. Quality guaranteed for 1 month if stored at 2 - 8°C.

Materials required but not provided

- Calibrated micropipettes and disposable tips.
- Distilled or de-ionized water.
- Plate washer: automatic or manual.
- Polypropylene tubes.
- Calibrated ELISA plate reader capable of measuring absorbance at 450 nm.
- Adhesive covers can be ordered separately. Please contact your local distributor.
- Centrifuge for 1 ml tubes.

6. WARNINGS AND PRECAUTIONS

- For research use only, not for diagnostic or therapeutic use.
- This kit should only be used by qualified laboratory staff.
- Do not under any circumstances add sodium azide as preservative to any of the components.
- Do not use kit components beyond the expiration date.
- Do not mix reagents from different kits and lots. The reagents have been standardized as a unit for a given lot. Use only the reagents supplied by the manufacturer.
- The assay has been optimized for the indicated standard range. Do not change the standard range.
- Open vials carefully: vials are under vacuum.
- It is advised to spin down streptavidin-peroxidase tubes before use.
- Do not ingest any of the kit components.
- Kit reagents contain 2-chloroacetamide as a preservative. 2-Chloroacetamide is harmful in contact with skin and toxic if swallowed. In case of accident or if you feel unwell, seek medical advise immediately.
- The TMB substrate is light sensitive, keep away from bright light. The solution should be colorless until use.
- The stop solution contains 2% oxalic acid and can cause irritation or burns to respiratory system, skin and eyes. Direct contact with skin and eyes should be strictly avoided. If contact occurs, rinse immediately with plenty of water and seek medical advise.
- Incubation times, incubation temperature and pipetting volumes other than those specified may give erroneous results.
- Do not reuse microwells or pour reagents back into their bottles once dispensed.
- Handle all biological samples as potentially hazardous and capable of transmitting diseases.
- Hemolyzed, hyperlipemic, heat-treated or contaminated samples may give erroneous results.
- Use polypropylene tubes for preparation of standard and samples. Do not use polystyrene tubes or sample plates.
- The standard is of human origin. It has been tested for various viruses and found negative. Since no test method can offer complete assurance that infectious agents are absent, this reagent should be handled as any potentially infectious human serum or blood specimen. Handle all materials in contact with this reagent according to guide-lines for prevention of transmission of blood-borne infections.

7. SAMPLE PREPERATION

Collection and handling

Serum or plasma

Collect blood using normal aseptic techniques. Blood samples should be kept on ice. If serum is used, separate serum from blood after clotting at room temperature within 1 hour by centrifugation (1500xg at 4°C for 15 min). Transfer the serum to a fresh polypropylene tube.

If plasma is used, separate plasma from blood within 20 minutes after blood sampling by centrifugation (1500xg at 4°C for 15 min). Transfer the plasma to a fresh polypropylene tube.

Bronchoalveolar lavage fluid (BALF)

Perform BAL during bronchoscopy by standardized washing of the segment. Aspirate each aliquot of physiologic fluid immediately after inspiration. Discard the first aliquot of recovered BALF. Collect the BALF in polypropylene tubes and keep it on ice. Separate cells from BALF by centrifugation (500xg at 4°C for 5 min). Filter cell free BALF through a layer of gauze to remove mucus strands.

Urine

Collect urine using normal aseptic techniques. Centrifuge the urine to remove debris (1500xg at 4 °C for 15 min). Transfer urine to a fresh polypropylene tube.

Storage

Store samples below -20°C, preferably at -70°C in polypropylene tubes. Storage at -20°C can affect recovery of human Properdin. Use samples within 24 hours after thawing. Avoid multiple freeze-thaw cycles which may cause loss of human Properdin activity and give erroneous results.

Do not use hemolyzed, hyperlipemic, heat-treated or contaminated samples.

Before performing the assay, samples should be brought to room temperature (18 – 25°C) and mixed gently. Prepare all samples (controls and test samples) prior to starting the assay procedure. Avoid foaming.

Dilution procedures

Serum or plasma samples

Human Properdin can be measured accurately if serum or plasma samples are diluted at least 5000x with supplied dilution buffer in polypropylene tubes.

BALF samples

Human Properdin can be measured accurately if BALF samples are diluted at least 10x with supplied dilution buffer in polypropylene tubes.

Urine samples

Human Properdin can be measured accurately if urine samples are diluted at least 10x with supplied dilution buffer in polypropylene tubes.

CSF samples

Human Properdin can be measured accurately if CSF samples are diluted at least 5x with supplied dilution buffer in polypropylene tubes.

Comment regarding recommended sample dilution

The mentioned dilution for samples is a minimum dilution and should be used as a guideline. The recovery of human Properdin from an undiluted sample is not 100% and may vary from

sample to sample. When testing less diluted samples it is advisable to run recovery experiments to determine the influence of the matrix on the detection of human Properdin.

Do not use polystyrene tubes or sample plates for preparation or dilution of the samples.

Guideline for dilution of samples

Please see the table below for recommended sample dilutions. Volumes are based on a total volume of at least 230 µl of diluted sample, which is sufficient for one sample in duplicate in the ELISA. For dilution of samples we recommend to use at least 10 µl of sample.

	Dilution	Pre-dilution	Amount of sample or pre-dilution required	Amount of Dilution buffer required
1.	10x	Not necessary	25 µl (sample)	225 µl
2.	20x	Not necessary	15 µl (sample)	285 µl
3.	50x	Not necessary	10 µl (sample)	490 µl
4.	100x	Not necessary	10 µl (sample)	990 µl
5.	500x	Recommended: 10x (see nr.1)	10 µl (pre-dilution)	490 µl
6.	1000x	Recommended: 10x (see nr.1)	10 µl (pre-dilution)	990 µl
7.	2000x	Recommended: 20x (see nr.2)	10 µl (pre-dilution)	990 µl
8.	5000x	Recommended: 50x (see nr.3)	10 µl (pre-dilution)	990 µl

Table 3

8. REAGENT PREPARATION

Allow all the reagents to equilibrate to room temperature (20 – 25°C) prior to use. Return to proper storage conditions immediately after use.

Wash buffer

Prepare wash buffer by mixing 60 ml of 20x wash buffer with 1140 ml of distilled or de-ionized water, which is sufficient for 2 x 96 tests. In case less volume is required, prepare the desired volume of wash buffer by diluting 1 part of the 20x wash buffer with 19 parts of distilled or de-ionized water.

Dilution buffer

Prepare dilution buffer by mixing 20 ml of the 10x dilution buffer with 180 ml of distilled or de-ionized water, which is sufficient for 2 x 96 tests. In case less volume is required, prepare the desired volume of dilution buffer by diluting 1 part of the 10x dilution buffer with 9 parts of distilled or de-ionized water. Concentrated dilution buffer may contain crystals. In case the crystals do not disappear at room temperature within 1 hour, concentrated dilution buffer can be warmed up to 37°C. Do not shake the solution.

Standard solution

The standard is reconstituted by pipetting the amount of dilution buffer mentioned on the CoA in the standard vial. Use the standard vial as Tube 1 in Figure 1. Prepare each human Properdin standard in polypropylene tubes by serial dilution of the reconstituted standard with dilution buffer as shown in Figure 1*. After reconstitution the standard cannot be stored for repeated use.

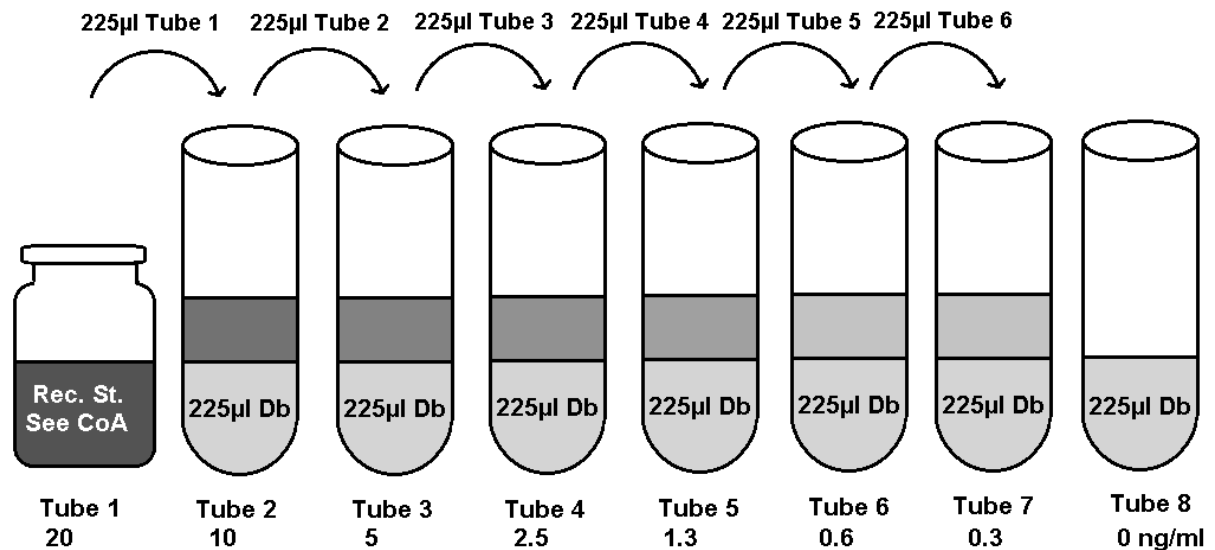


Figure 1

*) CoA: Certificate of Analysis, Rec. St: Reconstituted Standard, Db: Dilution buffer

Tracer solution

The tracer is reconstituted by pipetting 1 ml distilled or de-ionized water. Dilute the reconstituted 1 ml tracer with 11 ml dilution buffer, which is sufficient for 1 x 96 tests. In case less volume is required, prepare the desired volume of tracer by diluting 1 part of the reconstituted tracer with 11 parts of dilution buffer.

Streptavidin-peroxidase solution

It is advised to spin down streptavidin-peroxidase tubes before use. Prepare the streptavidin-peroxidase solution by mixing 0.25 ml of the 100x streptavidin-peroxidase solution with 24.75 ml dilution buffer, which is sufficient for 2 x 96 tests. In case less volume is required, prepare the desired volume of streptavidin-peroxidase solution by diluting 1 part of the 100x streptavidin-peroxidase solution with 99 parts of dilution buffer.

9. ELISA PROTOCOL

Bring all reagents to room temperature (20 - 25°C) before use.

1. Determine the number of test wells required, put the necessary microwell strips into the supplied frame, and fill out the data collection sheet. Return the unused strips to the storage bag with desiccant, seal and store at 2 - 8°C.
2. Transfer 100 µl in duplicate of standard, samples, or controls into appropriate wells. Do not touch the side or bottom of the wells.
3. Cover the tray and tap the tray to eliminate any air bubbles. Be careful not to splash liquid onto the cover.
4. Incubate the strips or plate for 1 hour at room temperature.
5. Wash the plates 4 times with wash buffer using a plate washer or as follows*:
 - a. Carefully remove the cover, avoid splashing.
 - b. Empty the plate by inverting plate and shaking contents out over the sink, keep inverted and tap dry on a thick layer of tissues.
 - c. Add 200 µl of wash buffer to each well, wait 20 seconds, empty the plate as described in 5b.
 - d. Repeat the washing procedure 5b/5c three times.
 - e. Empty the plate and gently tap on thick layer of tissues.
6. Add 100 µl of diluted tracer to each well using the same pipetting order as applied in step 2. Do not touch the side or bottom of the wells.
7. Cover the tray and incubate the tray for 1 hour at room temperature.
8. Repeat the wash procedure described in step 5a-e.
9. Add 100 µl of diluted streptavidin-peroxidase to each well, using the same pipetting order as applied in step 2. Do not touch the side or bottom of the wells.
10. Cover the tray and incubate the tray for 1 hour at room temperature.
11. Repeat the wash procedure described in step 5a-e.
12. Add 100 µl of TMB substrate to each well, using the same pipetting order as applied in step 2. Do not touch the side or bottom of the wells.
13. Cover the tray and incubate the tray for 30 minutes at room temperature. It is advised to control the reaction on the plate regularly. In case of strong development the TMB reaction can be stopped sooner. Avoid exposing the microwell strips to direct sunlight. Covering the plate with aluminium foil is recommended.
14. Stop the reaction by adding 100 µl of stop solution with the same sequence and timing as used in step 12. Mix solutions in the wells thoroughly by gently swirling the plate. Gently tap the tray to eliminate any air bubbles trapped in the wells.
15. Read the plate within 30 minutes after addition of stop solution at 450 nm using a plate reader, following the instructions provided by the instrument's manufacturer.

*) In case plate washer is used, please note: use of a plate washer can result in higher background and decrease in sensitivity. We advise validation of the plate washer with the manual procedure. Make sure the plate washer is used as specified for the manual method.

10. INTERPRETATION OF RESULTS

- Calculate the mean absorbance for each set of duplicate standards, control and samples.
- If individual absorbance values differ by more than 15% from the corresponding mean value, the result is considered suspect and the sample should be retested.
- The mean absorbance of the zero standard should be less than 0.3.
- Create a standard curve using computer software capable of generating a good curve fit. The mean absorbance for each standard concentration is plotted on the vertical (Y) axis versus the corresponding concentration on the horizontal (X) axis (logarithmic scale).
- If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.
- Samples that give a mean absorbance above the absorbance for the highest standard concentration are out of range of the assay. These samples should be retested at a higher dilution.

11. TECHNICAL HINTS

- User should be trained and familiar with ELISA assays and test procedure.
- If you are not familiar with the ELISA technique it is recommended to perform a pilot assay prior to evaluation of your samples. Perform the assay with a standard curve only following the instructions.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Completely empty wells before dispensing wash buffer, fill with wash buffer as indicated for each cycle and do not allow wells to sit uncovered or dry for extended periods.
- Since exact conditions may vary from assay to assay, a standard curve must be established for every run. Samples should be referred to the standard curve prepared on the same plate.
- Do not mix reagents from different batches, or other reagents and strips. Remainders should not be mixed with contents of freshly opened vials.
- Each time the kit is used, fresh dilutions of standard, sample, tracer, streptavidin-peroxidase and buffers should be made.
- Caps and vials are not interchangeable. Caps should be replaced on the corresponding vials.
- To avoid cross-contaminations, change pipette tips between reagent additions of each standard, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
- Waste disposal should be performed according to your laboratory regulations.

Technical support

Do not hesitate to contact our technical support team at support@hycultbiotech.com for inquiries and technical support regarding the human Properdin ELISA.

Hycult Biotech, Frontstraat 2a, 5405 PB Uden, the Netherlands

T: +31 (0)413 251 335, F: +31 (0)413 248 353

12. QUALITY CONTROL

The Certificate of Analysis included in this kit is lot specific and is to be used to verify results obtained by your laboratory. The absorption values provided on the certificate of analysis are to be used as a guideline only. The results obtained by your laboratory may differ.


This assay is designed to eliminate interference by soluble receptors, binding proteins, and other factors present in biological samples. Until all factors have been tested in the immunoassay, the possibility of interference cannot be excluded.

For optimal performance of this kit, it is advised to work according to good laboratory practice.

13. PERFORMANCE CHARACTERISTICS

Linearity

The linearity of the assay was determined by serially diluting a sample containing human Properdin. The diluted samples were measured in the assay.

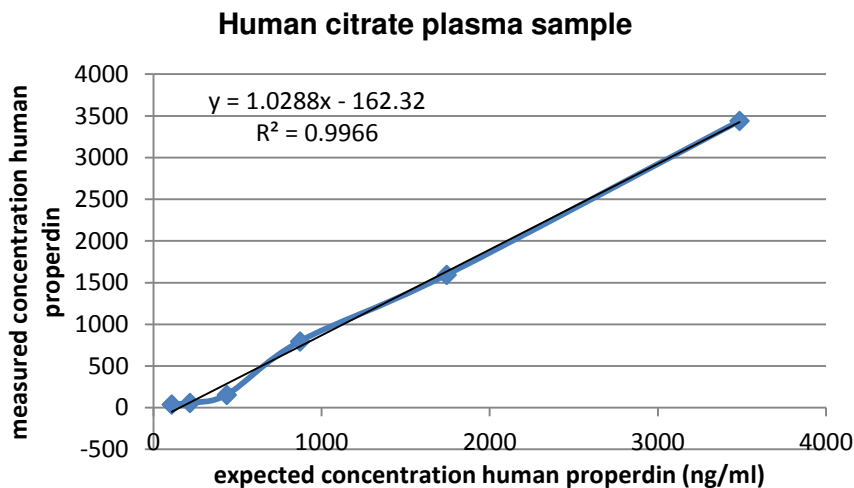

Recovery

Normal human blood samples (plasma) containing baseline levels of human Properdin, were spiked with human Properdin, in concentrations of 500 and 2000 ng/ml. Samples with and without human Properdin, were incubated for 1 hour at room temperature. Samples were measured using the ELISA. Values for human Properdin, ranged between 85.1% and 115.3%.

14. TROUBLESHOOTING

Warranty claims and complaints in respect of deficiencies must be logged before expiry date of the product. A written complaint containing lot number of the product and experimental data shall be sent to support@hycultbiotech.com.

Suggestions summarized below in Table 4 can be used as guideline in case of unexpected assay results.

Low absorbance	High absorbance	Poor duplicates	All wells positive	All wells negative	Possible cause
•	•		•	•	Kit materials or reagents are contaminated or expired
•					Incorrect reagents used
•		•	•		Lyophilized reagents are not properly reconstituted
•	•	•	•	•	Incorrect dilutions or pipetting errors
•		•			Improper plastics used for preparation of standard and/or samples
•	•				Improper incubation times or temperature
		•			Especially in case of 37°C incubation: plates are not incubated uniformly
•					Assay performed before reagents had reached room temperature
•	•	•	•	•	Procedure not followed correctly
				•	Omission of a reagent or a step
		•			Poor mixing of samples
	•		•		Low purity of water
	•	•			Strips were kept dry for too long during/after washing
	•	•	•		Inefficient washing
	•	•			Cross-contamination from other samples or positive control
		•	•		TMB solution is not clear or colourless
•	•				Wrong filter in the microtiter reader
	•	•			Airbubbles
		•			Imprecise sealing of the plate after use
•					Wrong storage conditions
•					Lamp in microplate reader is not functioning optimally

Table 4

15. REFERENCES

1. Harboe, M et al; The Role of Properdin in Zymosan- and Escherichia coli-Induced Complement Activation. *J Immunol* 2012, *189*: 2606